




Specifications

Model	B-V233P-A/B	B-V243Q-A/B	B-V433P-A/B	B-V443Q-A/B
Maximum balloon diameter	15 mm	15 mm	20 mm	20 mm
Sheath diameter	Distal end: 5.5 Fr		Proximal end: 7 Fr	
Sheath design	Over the wire	Distal wireguided	Over the wire	Distal wireguided
Working length	1900 mm	1950 mm	1900 mm	1950 mm
Minimum channel size	2.8 mm	3.2 mm	3.2 mm	3.2 mm
Compatible guidewire	0.89 mm (0.035 inch)			
Radiopaque band	One band at the distal end and one band at the proximal end of balloon			
Syringes	Three syringes are contained in the balloon package			
Lumen	8.5 mm/11.5 mm/15 mm	8.5 mm/11.5 mm/15 mm	15 mm/18 mm/20 mm	15 mm/18 mm/20 mm
Injection port	Triple lumen A: Above type B: Below type			

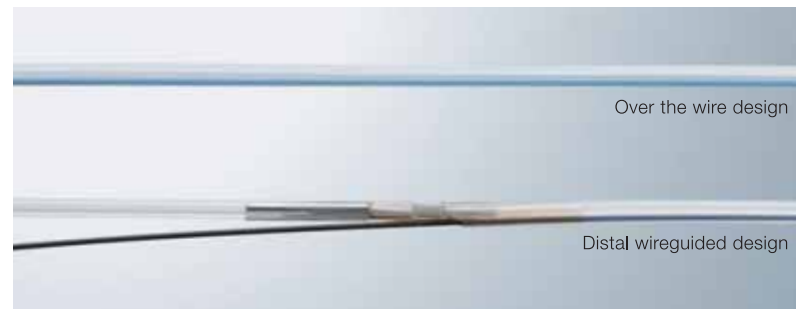

**OLYMPUS**

Your Vision, Our Future

OLYMPUS  
**EndoTherapy**  
BEYOND THE SCOPE

Single Use 3-Lumen Extraction Balloon V




Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

**OLYMPUS**

OLYMPUS MEDICAL SYSTEMS CORP.  
Shinjuku Monolith, 2-3-1 Nishi-Shinjuku, Shinjuku-ku, Tokyo 163-0914, Japan

For a complete listing of sales and distribution locations visit:  
[www.olympus.com](http://www.olympus.com)

Printed in Japan G000000-0000

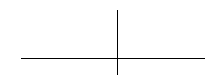

## Excellent balloon durability combined with efficient stone extraction capability provides choice and flexibility for ERCP physicians





### Excellent balloon durability

Thick, high-quality latex increases the durability and reduces the possibility of the balloon bursting during stone extraction.

Excellent balloon durability supports repeated stone extraction.



### Multiple balloon sizing

Balloons can be inflated to one of three diameters: 8.5 mm, 11.5 mm, and 15 mm or 15 mm, 18 mm, and 20 mm. Each balloon's size can be easily adjusted to suit the anatomical conditions of each case.

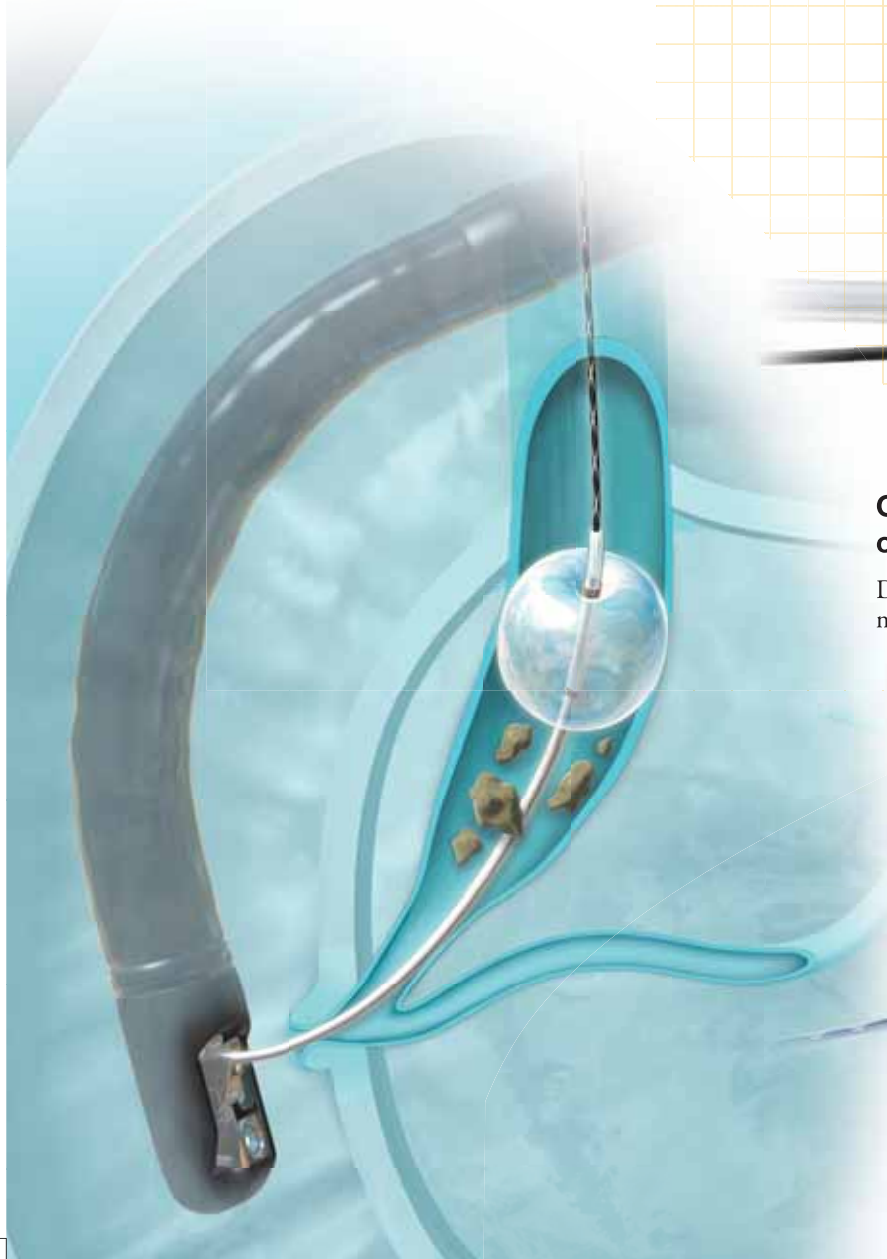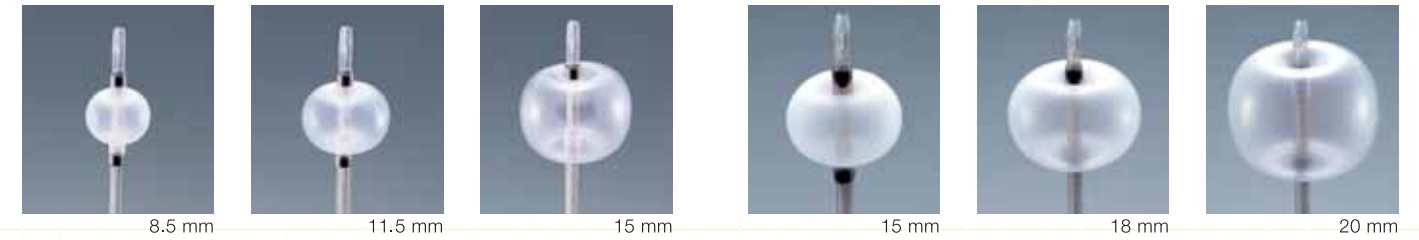

### Choice of distal wireguided or over the wire design

Distal wireguided option facilitates extraction of multiple stones with the guidewire still in the bile duct.

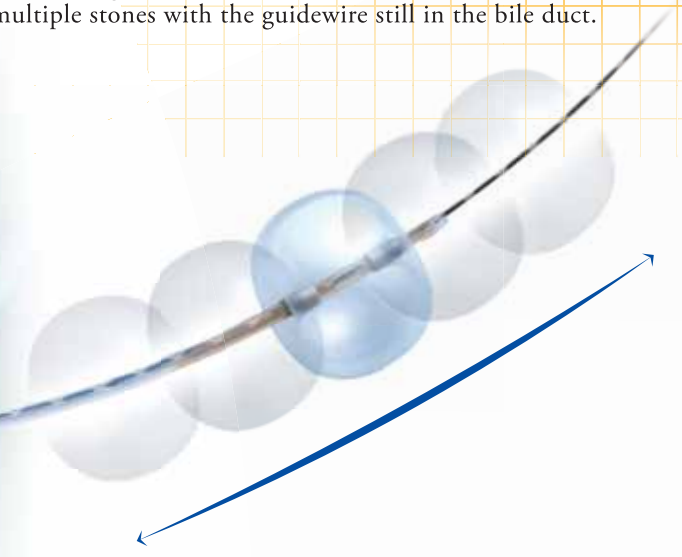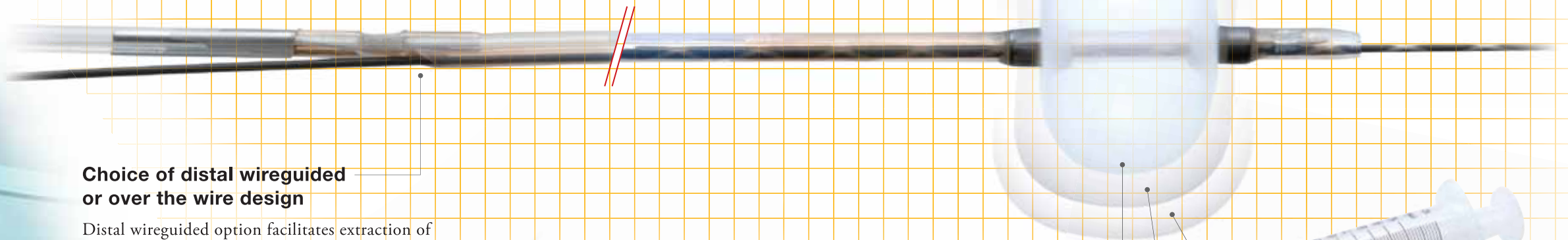

### Three pre-measured syringes for reliable inflation

Clearly marked syringes enable precise inflation of each balloon quickly without having to visually check the markings on the syringe.

