

Single Use Retrieval Basket V



Single Use Retrieval Basket V




# Versatile C-Hook Enhances Stone Management on Rotatable or Wire-guided Baskets

Offering unprecedented control, versatility, high performance, and a timesaving design that enables you to replace or switch devices in just seconds, the innovative V-System allows endoscopists to carry out a variety of ERCP procedures without assistance



## C-Hook

Attached to the scope control section, the convenient C-Hook makes it easy for you to maneuver the guidewire, inject contrast media, and manipulate the device handle — all while keeping a grip on the scope control section.

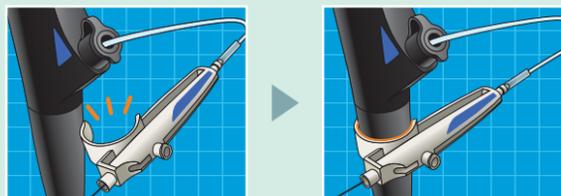

## V-Marking

Located on the proximal side of the sheath, the exclusive V-Marking lets you know when you can lower the V-Groove forceps elevator by showing you when the device tip has reached the distal end of the scope. Similarly, during withdrawal, it shows you when to raise the elevator to lock the guidewire.



## Unique Distal Eight-wire Basket is Ideal for Retrieval of Small Stones and with Four Proximal Wires to Facilitate Release of Captured Stones

### Eight-wire basket designed for the capture and removal of small stones.

The FlowerBasketV has eight wires on the distal side and four wires on the proximal side. The eight-wire section allows the precise retrieval of small stones while the four-wire section makes it easier to release them.



### Wire-guided ability for outstanding insertion capability

The FG-V431P and FG-V432P are designed to be wire guided with the wire passing only through the distal end of the basket. This makes it possible to insert the basket into a small incision on the papilla, and allows selective insertion into the intrahepatic bile duct. In addition, these baskets have the same expansion capability as our conventional models (FG-401Q/FG-22Q-1, etc.) because the guidewire is only attached to the basket's distal tip.



### Four-wire basket design for versatile retrieval

For versatile retrieval of a wide range of stones, the TetraCatchV features a standard four-wire retrieval basket configuration.

### Rotation mechanism and bullet-shaped tip

The rotatable basket (FG-V421PR/FG-V422PR) allows the retrieval of small stones, especially floating stones or stones located in the lower bile duct above the papilla. The bullet-shaped tip ensures smooth insertion into the bile duct.

Injection port for contrast media facilitates fluoroscopic visualization. Sterile single-use design for convenience both before and after procedures. Siliconized basket wire helps physicians to open and close the basket smoothly and facilitates "selective control" of the device with the C-Hook.

## Specifications

Model	Basket type	Rotatable/Wire-guided type	Working length	Basket opening width	Minimum channel size	Compatible guidewire
FG-V421PR	8 wire type	Rotatable type	1900 mm	20 mm	2.8 mm	0.89 mm (0.035 inch)
FG-V422PR	4 wire type			22 mm		
FG-V431P	8 wire type	Wire-guided type		20 mm	3.7 mm	
FG-V432P	4 wire type			22 mm		



### BML-110A-1

An emergency device, this is only to be used in the unlikely event of basket impaction due to an extremely hard calculus.

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

**OLYMPUS**<sup>®</sup>

www.olympus.com

**OLYMPUS MEDICAL SYSTEMS CORP.**  
Shinjuku Monolith, 3-1 Nishi-Shinjuku 2-chome, Shinjuku-ku, Tokyo 163-0914, Japan  
**OLYMPUS MEDICAL SYSTEMS EUROPA GMBH**  
Postfach 10 49 08, 20034 Hamburg / Wendenstrasse 14-18, 20097 Hamburg, Germany  
**OLYMPUS AMERICA INC.**  
3500 Corporate Parkway, P.O. Box 610, Center Valley, PA 18034-0610, U.S.A.  
**OLYMPUS LATIN AMERICA, INC.**  
5301 Blue Lagoon Drive, Suite 290 Miami, Florida 33126-2097, U.S.A.  
**KEYMED LTD.**  
Keymed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England

**OLYMPUS SINGAPORE PTE LTD.**  
491B River Valley Road #12-01/04, Valley Point Office Tower, Singapore 248373  
**OLYMPUS HONG KONG AND CHINA LIMITED**  
Room 1520-1527, Ocean Centre, 5 Canton Road, Tsimshatsui, Kowloon, Hong Kong  
**OLYMPUS (BEIJING) SALES & SERVICE CO., LTD.**  
Room.1202, NCI Tower, A12 Jianguomenwai Avenue Chaoyang District Beijing, 100022, China  
**OLYMPUS MOSCOW LIMITED LIABILITY COMPANY**  
117071, Moscow, Malaya Kaluzhskaya 19, bld. 1, fl.2, Russia  
**OLYMPUS AUSTRALIA PTY LTD**  
31 Gilby Road, Mount Waverley, VIC., 3149, Australia

Printed in Japan F1091SB-3-1007