


OLYMPUS[®]

Your Vision, Our Future

OLYMPUS
EndoTherapy
BEYOND THE SCOPE

EndoTherapy for small intestinal procedures

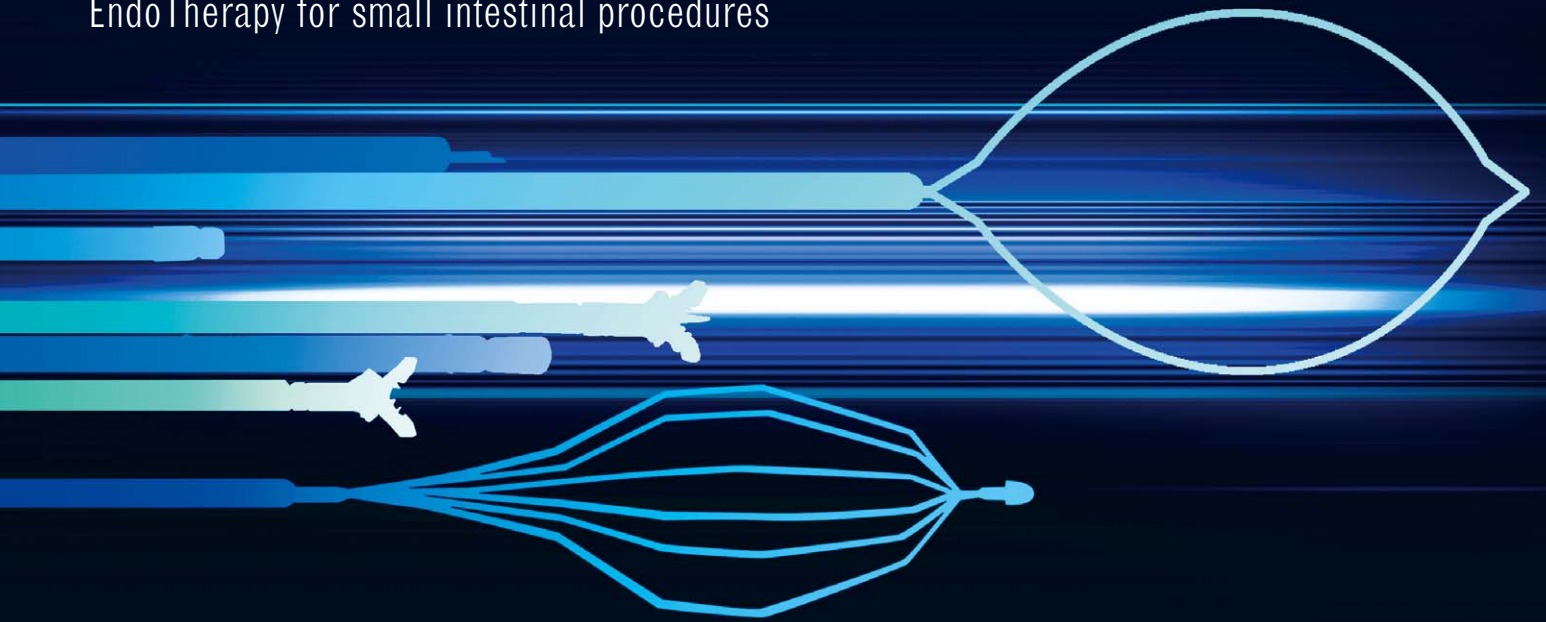

Versatile EndoTherapy accessories from Olympus offer exceptional support for small intestinal procedures, advancing enteroscopy into the future

Sampling

Biopsy Forceps



Single Use

Employment of an improved sheath section enables smooth insertion into a scope channel. These forceps are especially designed to prevent damage to the interior of a scope channel. Moreover, the newly developed blades (double-step blades) on the cups ensure sharp cutting quality, and swing jaw mechanism enables easier cup positioning for biopsy sites approached tangentially.



- FB-210U.A** Alligator Jaw without needle
- FB-220U.A** Alligator Jaw with needle
- FB-230U.A** Standard oval without needle
- FB-240U.A** Standard oval with needle

Minimum Channel Size	2.8 mm
Working Length	2300 mm
Quantity	20 pcs / Box

Reusable

A hole of approximately 1 mm diameter at the center of the cups designed to help minimize tissue damage and provide large sample yield.



- FB-28Y-1** Standard oval without needle
- FB-28Z-1** Standard oval without needle

Minimum Channel Size	2.8 mm
Working Length	FB-28Y-1 : 2700 mm FB-28Z-1 : 3000 mm
Quantity	1 pc / Box

Foreign Body Retrieval

Retrieval Basket

Especially designed for foreign body removal in small intestine, this retrieval basket offers reliable performance in retrieving resected tissue and other foreign objects with 6 wire basket and rotatable function.

Single Use



- FG-460YR** 6 wire basket

Minimum Channel Size	2.8 mm
Working Length	2700 mm
Basket Design	6 wire basket
Rotatability	Yes
Quantity	1 pc / Box

Polypectomy

Hot Biopsy Forceps



Disposable hot biopsy forceps with the slim outer diameter, the proprietary sheath design and improved flexibility ensure smooth insertion into the endoscope channel.

Single Use



- FD-210U** Alligator Jaw
- FD-230U** Standard oval

Minimum Channel Size	2.8 mm
Working Length	2300 mm
Quantity	5 pcs / Box

Electrosurgical Snare



Single Use

Disposable design with integrated handle for ease of use.



- SD-221U-25** Crescent type

Minimum Channel Size	2.0 mm
Working Length	2300 mm
Loop Opening Width	25 mm
Quantity	10 pcs / Box

Reusable

Especially designed for small intestinal procedure.



- SD-310Z-25** Oval type

Minimum Channel Size	2.8 mm
Working Length	3000 mm
Loop Opening Width	25 mm
Quantity	2 Wires and 2 tube sheaths




With the introduction of an innovative, yet surprisingly simple Single Balloon Enteroscope System, Olympus has literally opened up the small intestine to enteroscopy, helping this promising field take a quantum leap forward. Our sophisticated line of versatile EndoTherapy accessories combine with the advanced capabilities of the Single Balloon Enteroscope to create a synergy that will put you at the cutting edge of enteroscopy.

Hemostasis

Rotatable Clip Fixing Device **QuickClip2**

Especially designed for enteroscopy procedures, the enhanced rotation feature facilitates easy and effective targeting of lesions in the small intestines.

Single Use

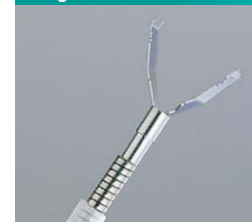

- HX-201YR-135**

Minimum Channel Size	2.8 mm
Working Length	2700 mm
Clip Jaw Angle	135°
Rotatability	Yes
Clip Opening Width	Standard
Quantity	5 pcs / Box

Bipolar Probes



Thanks to a specially formulated coating designed to minimize sticking to coagulated tissue and enhance hemostasis, these probes deliver bipolar hemostasis for uniform coagulation. A silver tip provides better conductivity than other metals.

Single Use



- CD-B420LA** Single plug
- CD-B422LA** Dual plug

Silver, Hydrophilic Tip Coating	
Minimum Channel Size	2.8 mm (CD-B420/B422-LA)
Working Length	3500 mm
Quantity	1 pc / Box

Injection Needle



Thanks to improved needle extrusion capability, enhanced click response, and excellent liquid feeding, the InjectorForce Max™ maximizes endoscopic injection performance.

Single Use



- NM-400Y-0423** Short bevel


Minimum Channel Size	2.8 mm
Working Length	2700 mm
Needle Length	4 mm
Gauge	23 G
Quantity	5 pcs / Box

Other Accessories

Washing Pipe / Spray Catheters

Especially designed for spraying dye, these catheters enable circular spraying of the liquid from the distal end.

Reusable

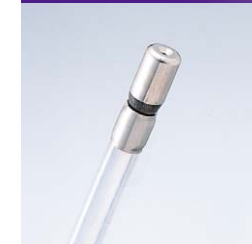

- PW-1V-1** Standard type
- PW-5V-1** Spray type


Minimum Channel Size	2.8 mm
Working Length	2400 mm
Feature	PW-1V-1 : Inside Chip (Routine washing) PW-5V-1 : Screw tip (Staining)
Quantity	PW-1V-1 : 2 pcs / Box (A), 1 pc / Box (B) PW-5V-1 : 1 pc / Box

Cleaning Brushes



Designed for one-time use, this brush is ideal for cleaning suction channels and openings of compatible Olympus endoscopes.


Single Use

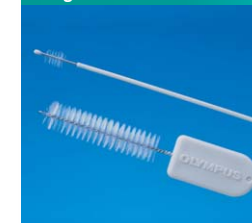

- BW-201T(A)** Single use channel cleaning brush
- BW-201T(B)** Channel cleaning brush and channel-opening cleaning brush

Working Length	2200 mm
Quantity	A : 50 pcs / Box B : 50 sets / Box

Reusable

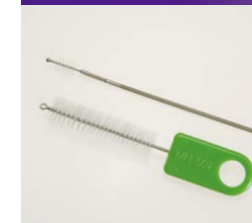


- BW-9Y** Coil sheath channel cleaning brush
- MH-507** Channel-port cleaning brush

Working Length	2850 mm (BW-9Y)
Quantity	1 pc / box

Distal Attachment

When attached to the distal end of a scope, this transparent cap assists in observation of the mucosa and helps achieve a good field of view.

Single Use

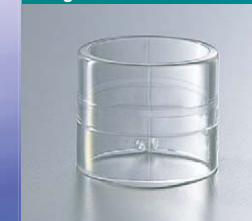

- D-201-10704** Soft Cap, Straight type

Max. Outer Diameter	11.0 mm
Distal Length from Scope	4 mm
Applicable Scope	SIF-Q180
Quantity	12 pcs / box

OLYMPUS

EndoTherapy
BEYOND THE SCOPE

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

OLYMPUS[®]

www.olympus.com

OLYMPUS MEDICAL SYSTEMS CORP.
Shinjuku Monolith, 3-1 Nishi-Shinjuku 2-chome, Shinjuku-ku, Tokyo 163-0914, Japan
OLYMPUS WINTER & IBE GMBH
Kuehnstrasse 61, 22045 Hamburg, Germany
OLYMPUS MEDICAL SYSTEMS EUROPA GMBH
Postfach 10 49 08, 20034 Hamburg / Wendenstrasse 14-18, 20097 Hamburg, Germany
OLYMPUS AMERICA INC.
3500 Corporate Parkway, P.O. Box 610, Center Valley, PA 18034-0610, U.S.A.
OLYMPUS LATIN AMERICA, INC.
5301 Blue Lagoon Drive, Suite 290 Miami, Florida 33126-2097, U.S.A.
OLYMPUS KEYMED GROUP LTD.
Keymed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England

OLYMPUS SINGAPORE PTE LTD.
491B River Valley Road #12-01/04, Valley Point Office Tower, Singapore 248373
OLYMPUS HONG KONG AND CHINA LIMITED
L43, Office Tower, Langham Place, 8 Argyle Street, Mongkok, Kowloon, Hong Kong
OLYMPUS (BEIJING) SALES & SERVICE CO., LTD.
Room 1202, NCI Tower, A12 Jianguomenwai Avenue Chaoyang District Beijing, 100022, China
OLYMPUS MOSCOW LIMITED LIABILITY COMPANY
117071, Moscow, Malaya Kaluzhskaya 19, bld. 1, fl.2, Russia
OLYMPUS AUSTRALIA PTY LTD
31 Gilby Road, Mount Waverley, VIC, 3149, Australia
OLYMPUS KOREA CO., LTD.
4F, 5F, Gyeongam Bldg., 157-27, Samsong-dong, Gangnam-ku, Seoul, Korea

F1227TT-0209