

OLYMPUS[®]




EVIS EXERA COLONOVIDEOSCOPE/SIGMOIDVIDEOSCOPE
OLYMPUS CF TYPE Q160L/I/S

TODAY'S MOST VERSATILE CHOICE FOR COLONOSCOPY

Offering sharp, clear images, large image size,
and auxiliary water capability for easy removal
of mucus and stool inside the colon,
this is today's most versatile scope
for colonoscopic examination and treatment




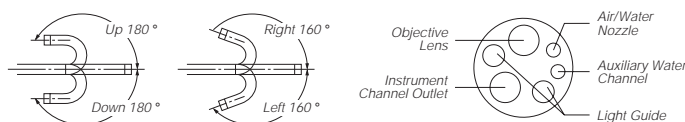
Enhanced Versatility and Outstanding Performance

- Superior image quality with large display size provides enhanced details throughout the image plane.
- When used in combination with our dedicated flushing pump and CV-160 video system center, a continuously clear view is ensured by an auxiliary water function that removes mucus and stool inside the colon at the touch of a switch on the scope. 
- Extra-wide 140° field of view and extensive 4-way angulation (180° up/down and 160° right/left).
- 3.7 mm diameter channel can handle a wide range of treatment requirements and provides high suction capability even when forceps are being used.
- Ergonomically designed grip enhances scope maneuverability while redesigned easy-to-access knobs and user-selectable switches improve operability.
- Fully compatible with the CV-160/140.
- Scope ID function stores individual scope information and displays it on the monitor to facilitate endoscopy suite management.




Specifications

Optical System	Field of view	140j
	Direction of view	0j Forward viewing
	Depth of field	3 to 100 mm
Distal End	Outer diameter	12.8 mm
Insertion Tube	Outer diameter	L/ I:12.8 mm, S:13.2 mm
Bending Section	Angulation range	Up 180j, Down 180j, Right 160j, Left 160j
Working Length		L:1680 mm, I:1330 mm, S:730 mm
Total Length		L:2000 mm, I:1650 mm, S:1040 mm
Instrument Channel	Inner diameter	3.7 mm
	Minimum visible distance	5 mm from the distal end
	Endo-Therapy accessory entrance/exit position in field of view	



Olympus business areas

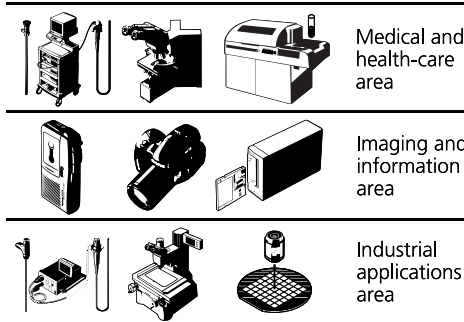

OLYMPUS
FOCUS ON LIFE

<http://www.olympus.com>

OLYMPUS OPTICAL CO., LTD.
San-Ei Building, 22-2, Nishi Shinjuku 1-chome, Shinjuku-ku, Tokyo, Japan
OLYMPUS OPTICAL CO. (EUROPA) GMBH.
Postfach 10 49 08, 20034 Hamburg / Wendenstrasse 14-16, 20097 Hamburg, Germany
OLYMPUS AMERICA INC.
2 Corporate Center Drive, Melville, N.Y. 11747-3157, U.S.A.
OLYMPUS KEYMED
Keymed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England
OLYMPUS SINGAPORE PTE LTD.
491B River Valley Road #12-01/04, Valley Point Office Tower, Singapore 248373
OLYMPUS BEIJING REPRESENTATIVE OFFICE
Rm. 818, South Tower, Beijing Kerry Centre, No.1 Guanghua Road, Chaoyang District, Beijing, 100020, China
OLYMPUS MOSCOW LIMITED LIABILITY COMPANY
117071, Moscow, Malaya Kaluzhskaya 19, bld. 1, fl.2, Russia
OLYMPUS AUSTRALIA PTY. LTD.
1/104 Ferntree Gully Road, Oakleigh, VIC 3166, Australia