

# OLYMPUS®



EVIS EXERA COLONOVIDEOSCOPE

## OLYMPUS CF TYPE 2T160L/I

# POWERFUL DUAL-CHANNEL PERFORMANCE

Featuring two independent channels with diameters of 3.7 mm and 3.2 mm, this scope offers the most powerful treatment and suction capability for colonoscopic procedures




## Main Features for Powerful Performance

- Dual channel design, with a 3.7 mm diameter channel for suction, and a 3.2 mm diameter channel for insertion of an Endo-Therapy accessory, ensures powerful suction capability and allows both channels to be used for Endo-Therapy accessories.
- 3 mm minimum depth of field for close-up observation.
- Extra-wide 140° field of view allows easier observation of a wider area of the site.
- 4-way angulation (180° up/down and 160° right/left) enables complete and comprehensive examination of the colon.
- 13.7 mm outer diameter insertion tube.
- Compatible with laser probes for expanded treatment versatility.
- Fully compatible with the CV-160/140/100.
- Scope ID function stores individual scope information and displays it on the monitor to facilitate endoscopy suite management.




## Specifications

Optical System	Field of view	140°
	Direction of view	0° Forward viewing
	Depth of field	3 to 100 mm
Distal End	Outer diameter	13.8 mm
	Insertion Tube	Outer diameter
Bending Section	Angulation range	Up 180°, Down 180°, Right 160°, Left 160°
Working Length		L:1680 mm, l:1330 mm
Total Length		L:2005 mm, l:1655 mm
Instrument Channel	Inner diameter	A:3.7 mm, B:3.2 mm
	Minimum visible distance	5 mm from the distal end
	Endo-Therapy accessory entrance/exit position in field of view	

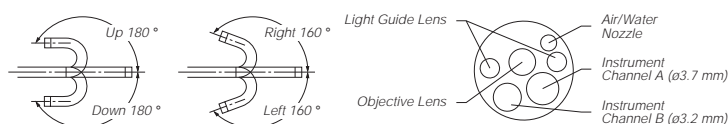

## Olympus business areas

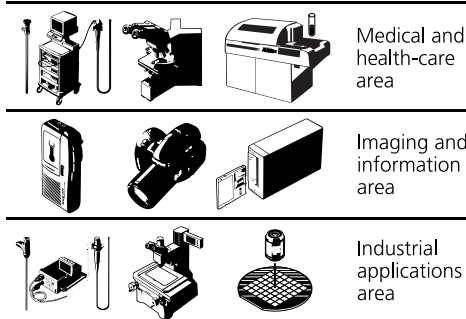

**OLYMPUS**  
FOCUS ON LIFE

**OLYMPUS OPTICAL CO., LTD.**  
2951 Ishikawa-cho, Hachioji-shi, Tokyo 192-8507, Japan  
**OLYMPUS OPTICAL CO. (EUROPA) GMBH.**  
Postfach 10 49 08, 20034 Hamburg / Wendenstrasse 14-18, 20097 Hamburg, Germany  
**OLYMPUS AMERICA INC.**  
2 Corporate Center Drive, Melville, N.Y. 11747-3157, U.S.A.  
**OLYMPUS KEYMED**  
Keymed House, Stovok Road, Southend-on-Sea, Essex SS2 5QH, England  
**OLYMPUS SINGAPORE PTE LTD.**  
491B River Valley Road #12-01/04 Valley Point Office Tower, Singapore 248373  
**OLYMPUS HONG KONG AND CHINA LIMITED.**  
Room 1520-1527, Ocean Centre, 5 Canton Road, Tsimshatsui, Kowloon, Hong Kong  
**OLYMPUS BEIJING REPRESENTATIVE OFFICE**  
Rm.818, South Tower, Beijing Kerry Centre, No.1 Guanghua Road, Chaoyang District, Beijing, 100020, China  
**OLYMPUS MOSCOW LIMITED LIABILITY COMPANY**  
117071, Moscow, Malaya Kaluzhskaya 19, bld. 1, fl.2, Russia  
**OLYMPUS AUSTRALIA PTY. LTD.**  
1/104 Ferntree Gully Road, Oakleigh, VIC 3166, Australia

<http://www.olympus.com>