

OLYMPUS[®]




EVIS EXERA GASTROINTESTINAL VIDEOSCOPE
OLYMPUS GIF TYPE 2T160

EXCEPTIONAL DUAL-CHANNEL VERSATILITY

Dual-channel design not only gives you more versatile treatment capability in the upper digestive tract, it also increases efficiency with enhanced suction capability




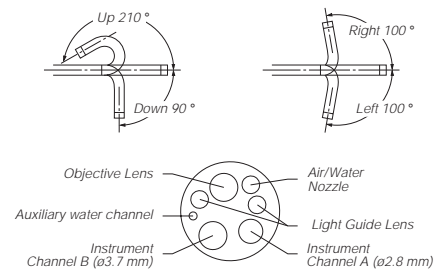
Main Features for Exceptional Performance

- Dual channel design, with a 3.7 mm diameter channel for suction, and a 2.8 mm diameter channel for insertion of an Endo-Therapy accessory, ensures powerful suction capability and allows, both channels to be used for Endo-Therapy accessories.
- When used in combination with our dedicated flushing pump and CV-160 video system center, an auxiliary water function is available to provide a continuously clear view at the touch of a switch on the scope. 
- 3 mm minimum depth of field for close-up observation.
- 12.6 mm outer diameter insertion tube.
- 120° field of view and 4-way angulation (210° up, 90° down and 100° right/left) enable comprehensive examination of the upper digestive tract.
- Compatible with laser probes for expanded treatment versatility.
- Fully compatible with the CV-160/140/100.
- Scope ID function stores individual scope information and displays it on the monitor to facilitate endoscopy suite management.




Specifications

Optical System	Field of view	120°
	Direction of view	0° Forward viewing
	Depth of field	3 to 100 mm
Distal End	Outer diameter	13.2mm
Insertion Tube	Outer diameter	12.6 mm
Bending Section	Angulation range	Up 210°, Down 90°, Right 100°, Left 100°
Working Length		1030 mm
Total Length		1350 mm
Instrument Channel	Inner diameter	A:2.8 mm, B:3.7 mm
	Minimum visible distance	4 mm from the distal end
	Endo-Therapy accessory entrance/exit position in field of view 	



Olympus business areas

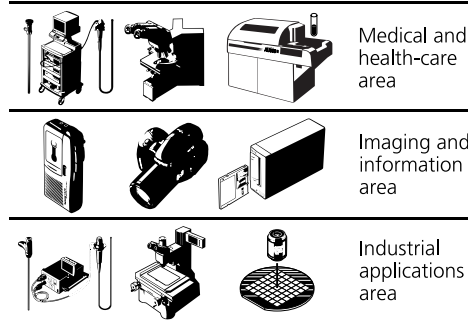

OLYMPUS
FOCUS ON LIFE

OLYMPUS OPTICAL CO., LTD.
2951 Ishikawa-cho, Hachioji-shi, Tokyo 192-8507, Japan
OLYMPUS OPTICAL CO. (EUROPA) GMBH.
Postfach 10 49 08, 20034 Hamburg / Wendenstrasse 14-18, 20097 Hamburg, Germany
OLYMPUS AMERICA INC.
2 Corporate Center Drive, Melville, N.Y. 11747-3157, U.S.A.
OLYMPUS KEYMED
Keymed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England
OLYMPUS SINGAPORE PTE LTD.
491B River Valley Road #12-01/04 Valley Point Office Tower, Singapore 248373
OLYMPUS HONG KONG AND CHINA LIMITED.
Room 1520-1527, Ocean Centre, 5 Canton Road, Tsimshatsui, Kowloon, Hong Kong
OLYMPUS BEIJING REPRESENTATIVE OFFICE
Rm.818, South Tower, Beijing Kerry Centre, No.1 Guanghua Road, Chaoyang District, Beijing, 100020, China
OLYMPUS MOSCOW LIMITED LIABILITY COMPANY
117071, Moscow, Malaya Kaluzhskaya 19, bld. 1, fl. 2, Russia
OLYMPUS AUSTRALIA PTY. LTD.
1/104 Ferntree Gully Road, Oakleigh, VIC 3166, Australia

<http://www.olympus.com>