


# OLYMPUS®

Your Vision, Our Future

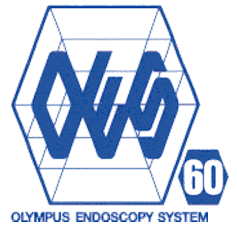

## OES Bronchofiberscope **OLYMPUS BF TYPE MP60**

**Slim design with a wide channel extends your reach  
into the bronchi**




With an outer diameter of only 4.0 mm at the distal end, this ingenious design incorporates a wide 2.0 mm channel



- Despite an outer diameter of only 4.0 mm at the distal end, this scope incorporates a large 2.0 mm instrument channel that can accommodate standard forceps for use in peripheral lung biopsies.
- Wide 120° field of view enables accurate observation of a wider area.
- 180° up / 130° down angulation for easier maneuvering in the bronchi.
- 600 mm working length ideal for use in the tracheobronchial tree.
- Compatible with both OES and EVIS video systems.

■ Specifications

Optical System	Field of view	120°
	Direction of view	0° (Forward viewing)
	Depth of field	3~50mm
Insertion Tube	Distal end outer diameter	4.4mm
	Insertion tube outer diameter	4.4mm
	Working length	600mm
Instrument Channel	Channel inner diameter	2.0mm
	Minimum visible distance	3mm from distal end
Bending Section	Angulation range	Up180°, Down130°

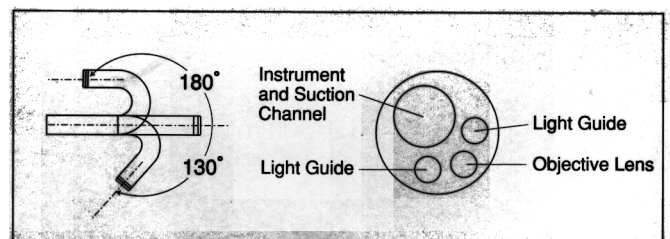

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

**OLYMPUS®**

OLYMPUS CORPORATION  
Shinjuku Monolith, 3-1 Nishi-Shinjuku 2-chome, Shjngjuku-ku, Tokyo 163-0914, Japan  
OLYMPUS WINTER & IBE GMBH  
Kaiserkrasse 61, 22045 Hamburg, Germany  
OLYMPUS EUROPA GMBH  
Postfach 718-49 08, 20334 Hamburg / Wendenstrasse 14-18, 20007 Hamburg, Germany  
OLYMPUS AMERICA INC.  
2 Corporate Center Drive, Marietta, N.Y. 11747-3157, U.S.A.  
OLYMPUS LATIN AMERICA, INC.  
3100 Blue Lagoon Drive, Suite 390 Miami, Florida 33126-2067, U.S.A.  
KEYMED LTD.  
Strydom House, Brook Road, Southend-on-Sea, Essex S82 5QH, England

OLYMPUS SINGAPORE PTE LTD.  
81B River Valley Road #12-01/04 Valley Point Office Tower, Singapore 246373  
OLYMPUS HONG KONG AND CHINA LIMITED  
Room 1505-1527, Ocean Centre, 6 Canton Road, Tsimshabau, Kowloon, Hong Kong  
OLYMPUS BEIJING REPRESENTATIVE OFFICE  
Room 1401, East Tower, Gong Yuan Hotel Royal Palace, Nuo Guo YuanXie, JianGuoMenwai, DongZheng District, Beijing, 100006, China  
OLYMPUS MOSCOW LIMITED LIABILITY COMPANY  
117071, Moscow, Malaya Prolskaya 19, Bld. 1, R.S., Russia  
OLYMPUS AUSTRALIA PTY LTD.  
11 Gilby Road, Mount Waverley, VIC. 3149, Australia